

2810 Clark Avenue • St. Louis, MO 63103-2574 • (314) 531-8080 • FAX (314) 531-8085
 Chemical, Metallurgical, Mechanical, Nondestructive, Environmental Testing, Analyses and Field Service.

RC IMAGING
 50 Old Hojack Lane
 Hilton, NY 14468



November 10, 2020
 Lab No. 20P-5262
 P.O. No. 1062601
 Page 1 of 6

Attention: Robert Wright

REPORT OF MECHANICAL TESTS

SAMPLE ID: RC Imaging PN 20-2407-700 U-Arm Two Step Platform
SUBJECT: Cyclic Load of HDPE Step
METHOD: Per Client Provided Load Test Plan and Drawing
TEST INSTRUMENT: Instron Model 5500R w/ Bluehill 3 Software
RESULTS: Specimen was not damaged during testing.

Sample ID	Cyclic Load For 250 Cycles (*Loading Plate is 40lbs) (lbs)	Final Displacement After Testing (in)
HDPE Step	710*	0.057

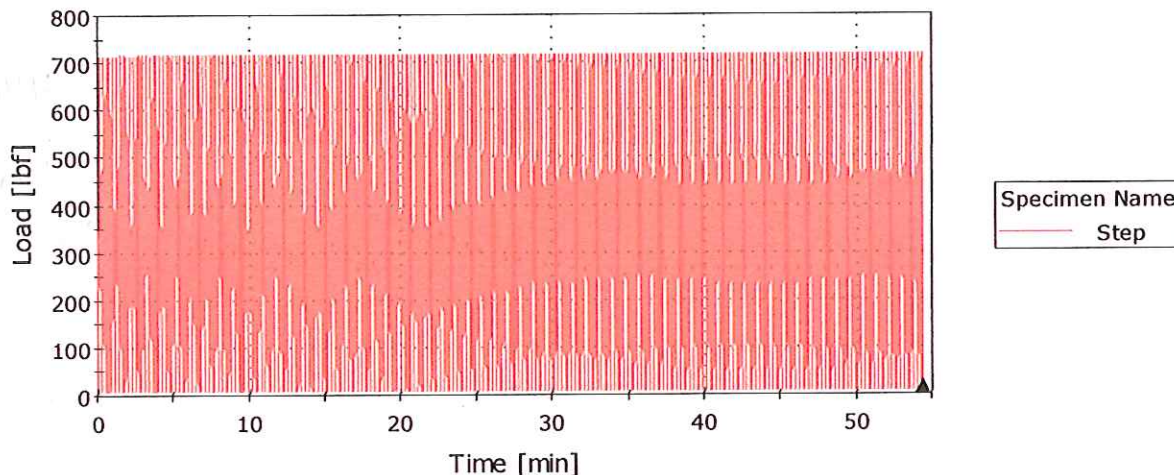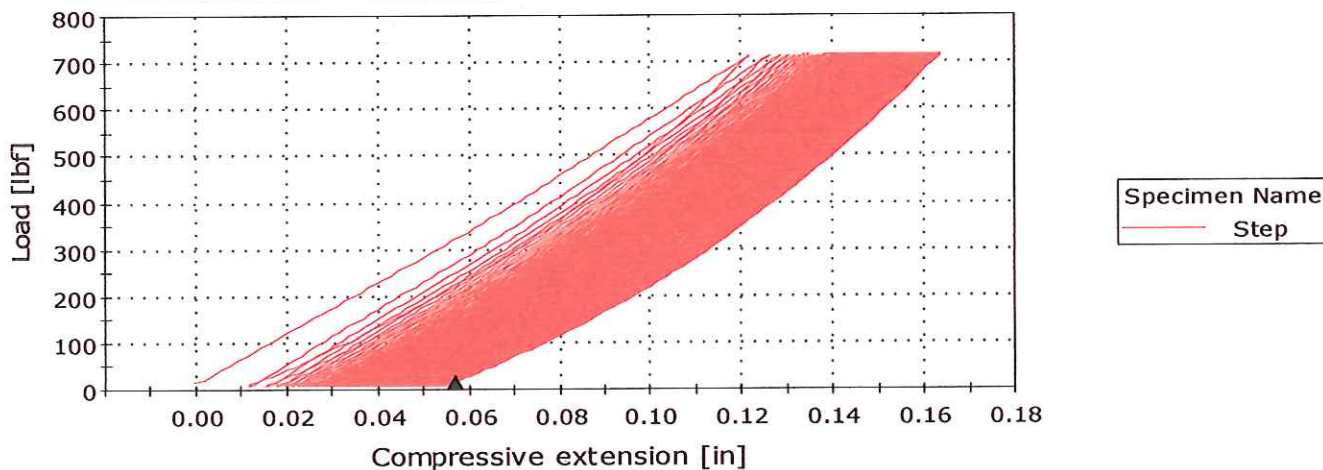

2810 Clark Avenue • St. Louis, MO 63103-2574 • (314) 531-8080 • FAX (314) 531-8085
 Chemical, Metallurgical, Mechanical, Nondestructive, Environmental Testing, Analyses and Field Service.

RC IMAGING
 50 Old Hojack Lane
 Hilton, NY 14468



November 10, 2020
 Lab No. 20P-5262
 P.O. No. 1062601
 Page 2 of 6

Attention: Robert Wright

REPORT OF MECHANICAL TESTS

SAMPLE ID: RC Imaging PN 20-2407-700 U-Arm Two Step Platform
SUBJECT: Static Load of HDPE Step
METHOD: Per Client Provided Load Test Plan and Drawing
TEST INSTRUMENT: Instron Model 5500R w/ Bluehill 3 Software
RESULTS: Specimen was not damaged during testing.

Sample ID	Static Load (*Loading Plate is 40lbs) (lbs)	Final Displacement After Testing (in)
HDPE Step	1460*	0.120

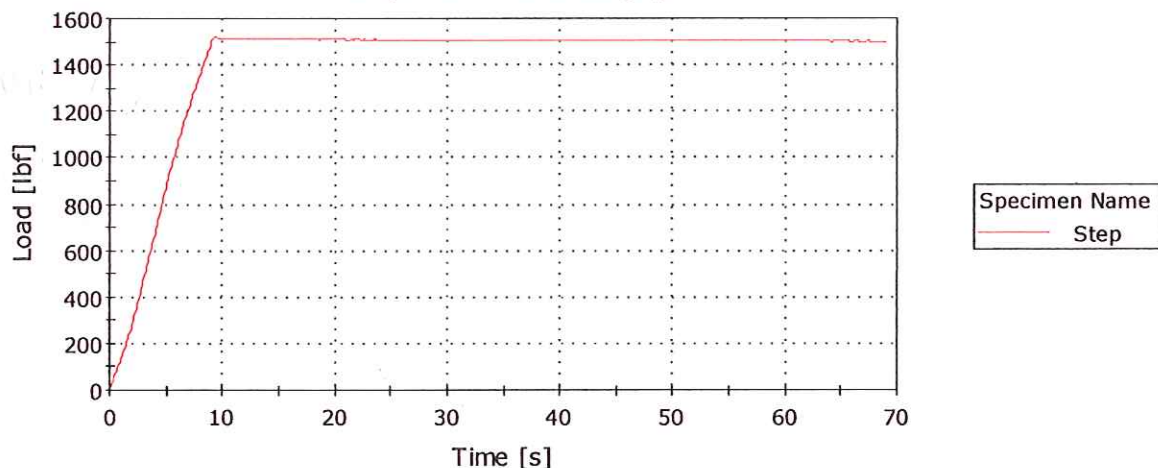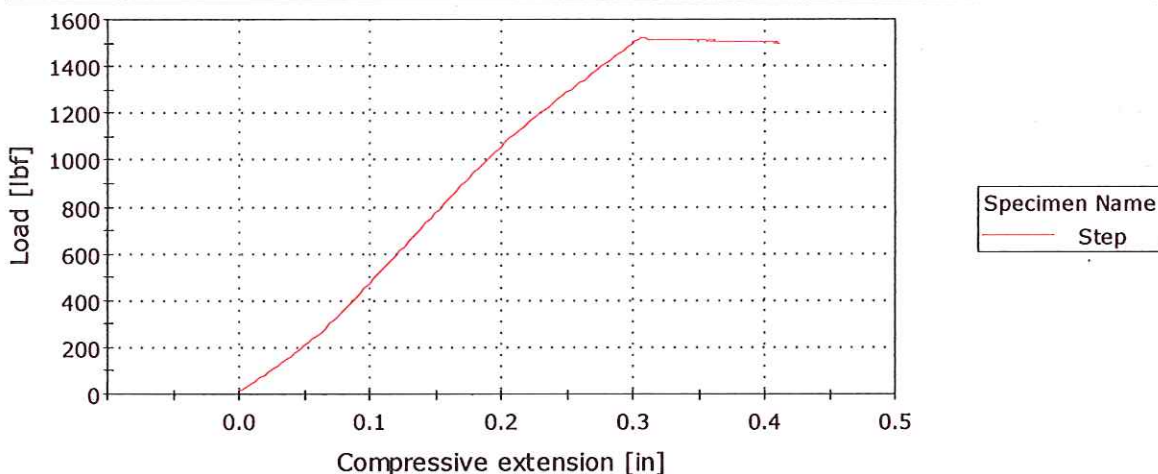

2810 Clark Avenue • St. Louis, MO 63103-2574 • (314) 531-8080 • FAX (314) 531-8085
Chemical, Metallurgical, Mechanical, Nondestructive, Environmental Testing, Analyses and Field Service.

RC IMAGING
50 Old Hojack Lane
Hilton, NY 14468

November 10, 2020
Lab No. 20P-5262
P.O. No. 1062601
Page 3 of 6

Attention: Robert Wright

REPORT OF MECHANICAL TESTS

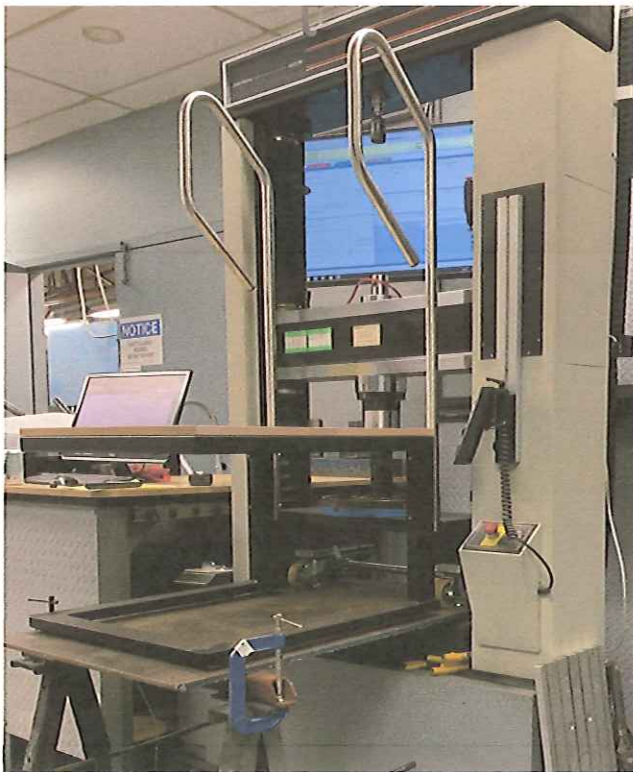


Figure 1: Test Setup For Step Load Tests



Figure 2: Test Setup For Step Load Tests

2810 Clark Avenue • St. Louis, MO 63103-2574 • (314) 531-8080 • FAX (314) 531-8085
 Chemical, Metallurgical, Mechanical, Nondestructive, Environmental Testing, Analyses and Field Service.

RC IMAGING
 50 Old Hojack Lane
 Hilton, NY 14468



November 10, 2020
 Lab No. 20P-5262
 P.O. No. 1062601
 Page 4 of 6

Attention: Robert Wright

REPORT OF MECHANICAL TESTS

SAMPLE ID: RC Imaging PN 20-2407-700 U-Arm Two Step Platform
SUBJECT: Cyclic Load of Acrylic Platform
METHOD: Per Client Provided Load Test Plan and Drawing
TEST INSTRUMENT: Instron Model 5500R w/ Bluehill 3 Software
RESULTS: Specimen was not damaged during testing.

Sample ID	Cyclic Load For 250 Cycles (*Loading Plate is 40lbs) (lbs)	Final Displacement After Testing (in)
Acrylic Platform	710*	0.011

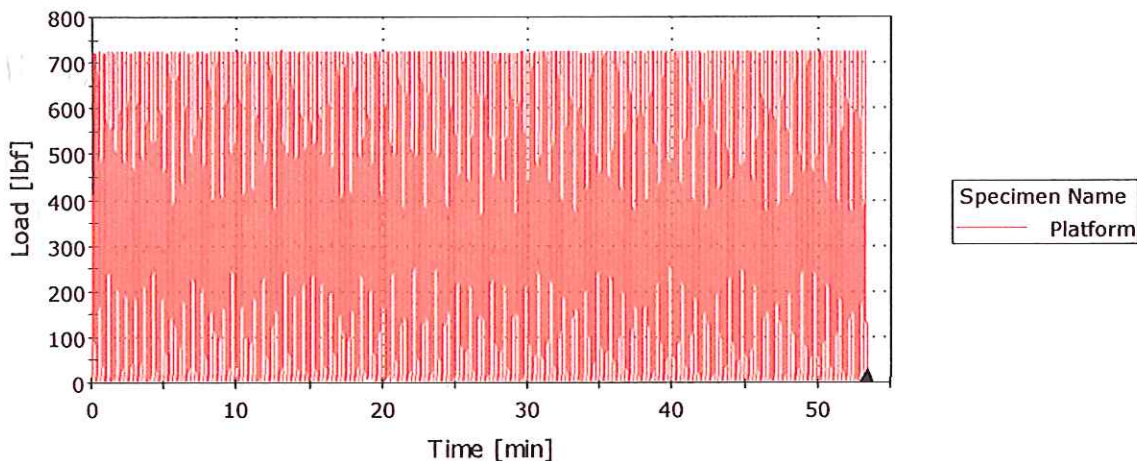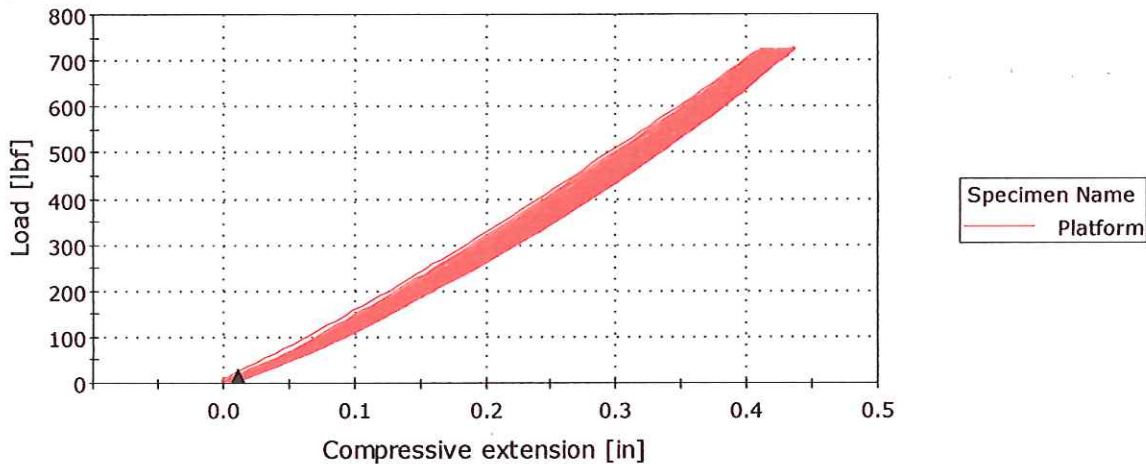

2810 Clark Avenue • St. Louis, MO 63103-2574 • (314) 531-8080 • FAX (314) 531-8085
 Chemical, Metallurgical, Mechanical, Nondestructive, Environmental Testing, Analyses and Field Service.

RC IMAGING
 50 Old Hojack Lane
 Hilton, NY 14468



November 10, 2020
 Lab No. 20P-5262
 P.O. No. 1062601
 Page 5 of 6

Attention: Robert Wright

REPORT OF MECHANICAL TESTS

SAMPLE ID: RC Imaging PN 20-2407-700 U-Arm Two Step Platform
SUBJECT: Static Load of Acrylic Platform
METHOD: Per Client Provided Load Test Plan and Drawing
TEST INSTRUMENT: Instron Model 5500R w/ Bluehill 3 Software
RESULTS: Specimen was not damaged during testing.

Sample ID	Static Load (*Loading Plate is 40lbs) (lbs)	Final Displacement After Testing (in)
Acrylic Platform	1460*	0.059

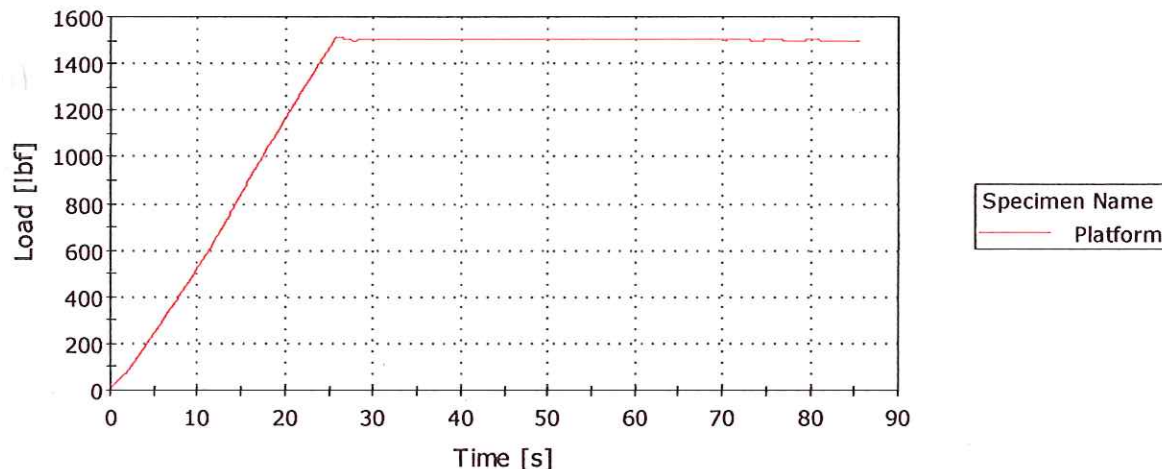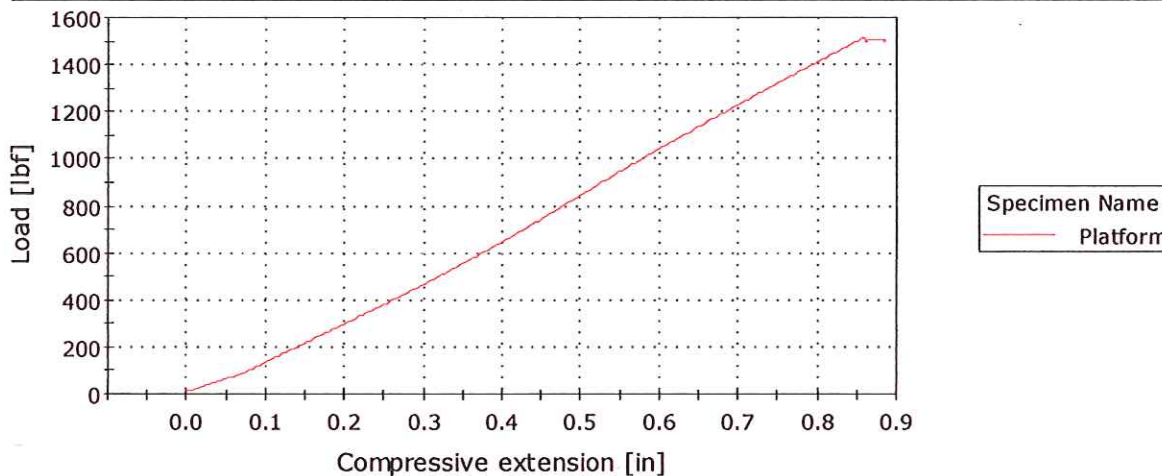

2810 Clark Avenue • St. Louis, MO 63103-2574 • (314) 531-8080 • FAX (314) 531-8085
Chemical, Metallurgical, Mechanical, Nondestructive, Environmental Testing, Analyses and Field Service.

RC IMAGING
50 Old Hojack Lane
Hilton, NY 14468

November 10, 2020
Lab No. 20P-5262
P.O. No. 1062601
Page 6 of 6

Attention: Robert Wright

REPORT OF MECHANICAL TESTS



Figure 3: Test Setup For Platform Load Tests Figure 4: Test Setup For Platform Load Tests

KS/ks


Karl Schmitz, Director
Materials Testing